


FUNCTIONAL NEUROSURGERY

# Single-use biopsy probes

With Sedan side-cutting window

**inomed** 

Closed window

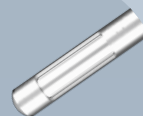

Open window



# Single-use biopsy probes


## With Sedan side-cutting window

### EFFICIENT BIOPSIES

The biopsy probes are designed to collect multiple tissue samples along a single trajectory from one or more target points. Each probe has a concentric inner and outer tube for precision and reliability.


#### Closed window

Ensures safe insertion and removal of the probe from brain tissue

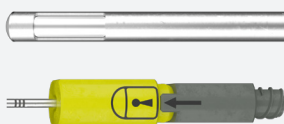


#### Open window

Enables targeted aspiration and collection of brain tissue

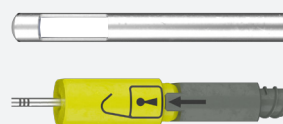

### PRECISE TARGETING

The target point is located at the centre of the side-cutting window, enabling accurate and controlled tissue sampling.

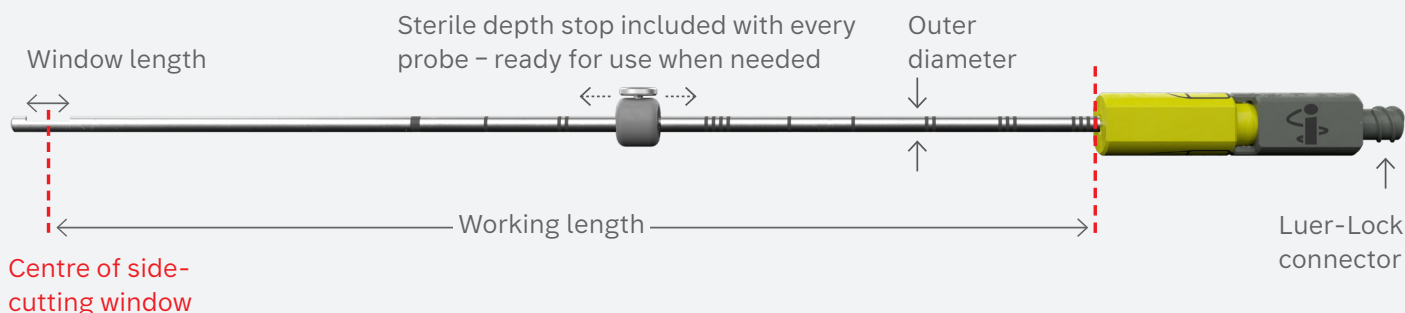

### FLEXIBILITY

inomed offers a broad selection of biopsy probes available in two working lengths, two window lengths and two outer diameters, providing versatile options for a variety of clinical requirements.

Article No.	Window length [mm]	Outer diameter [mm]	Colour
805610	5	2.1	Yellow
805611	8	2.1	Yellow
805612	5	2.5	Blue
805613	8	2.5	Blue

Article No.	Window length [mm]	Outer diameter [mm]	Colour
805620	5	2.1	Yellow
805621	8	2.1	Yellow
805622	5	2.5	Blue
805623	8	2.5	Blue